

HyBase 3000 Operating Table











Practical functionality, flexible versatility and high reliability

Functional Display of Remote Panel Longitudinal slide to head end Normal position Table down Trendelenburg position Reverse Trendelenburg position Back plate down Flex position Re-flex position



Build-in elevator for kidney bridge & chest bridge



Engineered for increasingly complex surgical scenes, HyBase needs least effort and time to position the patient in a safe, comfortable and efficient way. Insightful to the surgeon and nurse's needs, Mindray is answering demands for a total solution of operating table with:

- Better mattress for the comfort of patients
- Safe and efficient patient positioning
- Know-how for maximizing the versatility



- Multi-layer decompression pad with waterproof and seamless design
- Longitudinal slide up to 300 mm for free access to Carm
- Max. overall load up to 230kg in normal position
- Battery power for 50-80 operations
- Mechanical brake system for optimum stability
- Return to original position with one button
- Flex & Reflex position with one button
- Radiolucent table plate for free C-arm access
- Three control panels including remote panel, backup control & foot switch



Convenient positioning technique for all surgical disciplines

A successful operation starts with the excellent patient positioning, and thus requests convenient access, excellent results, time saving and efficiency. Outstanding flexibility of HyBase 3000 together with various accessories makes every possible patient positioning quickly and securely established.



Intervertebral operating position
Accessory:
0092-30-102639 Intervertebral disc operation pad
(Memorise foam, one piece)
115-005864-00 Prone position pillow
0092-30-102634 Arm board



Rectal operating position
Accessory:
115-007049-00 Rectal positioning device
115-005864-00 Prone position pillow
0092-30-102640 Body strap, one piece
0092-30-102634 Arm board



Arm or hand operating position Accessory: 115-007047-00 Arm operating table 115-005872-00 Head ring



General supine operating position with longitudinal slide of 300mm to head end
Accessory:
115-005872-00 Head ring

We use high quality stainless steel material in operating tables and accessories. In order to get the surgeons access to C arm, all the pads or plate are radiolucent for X-ray check. To avoid hurting the sensitive part of the patients, we provide Gel material pads for patients positioning which are similar as human tissue.



General thyroid operating position



Ophthalmic operating position
Accessory:
115-007045-00 Ophthalmic operating plate
115-005872-00 Head ring
0092-30-102642 Neurosurgical head support adapter



Neurosurgery operating position
Accessory:
115-012835-00 Adjustable horseshoe shaped head rest
0092-30-102642 Neurosurgical head support adapter



Neurosurgical adapter can connect Mayfield frame Accessory: 0092-30-102642 Neurosurgical head support adapter

Convenient positioning technique for all surgical disciplines

Developed in close collaboration with end users, the accessories fulfill all kind of surgery positioning demands and place the patients in an optimum way which brings benefits both for surgeons and patients.



Lithotomy operating position for gynecology & urology

115-007043-00 A pair of Leg holders with loops, with radial universal clamps

115-013018-00 Drainage bowl for Gynecology & Urology



Lithotomic operating position

0092-30-102635 one piece of heavy version leg support, with one piece of heavy radial universal clamp, pad & one piece of fasten belt

115-013018-00 Drainage bowl for Gynecology & Urology



Lateral position for MIS(minimally invasive surgery) MIS quickly est ablished with 120mm build-in elevator



Lateral position for MIS(minimally invasive surgery) Accessory:

0092-30-102634 Arm board

0092-30-102648 Lateral position arm support

0092-30-102638 Large body support

0092-30-102669 Small body support

115-005811-00 115-005872-00 Head ring

Arm/leg pad



Femoral operating position

Accessory:

115-009425-00 Orthopedic extension device with transport Cart



115-007057-00 Lateral femoral extension assembly 115-009425-00 orthopedic extension pre-sele cted



Accessory:

115-007055-00 Tibia extension assembly

115-009425-00 orthopedic extension pre-selected



Tibia operating position

Accessory: 115-007055-00 Tibia extension assembly 115-009425-00 Orthopedic extension device with transport Cart



Traumatology for the neck of femurs

Accessory:

115-007055-00 Tibia extension assembly 115-009425-00 Orthopedic extension device with transport Cart

All the advanced and reliable strengths from Mindray technology



With the new SFC (Special Foam Core) cushions, HyBase 3000 provides excellent pressure distribution, significantly avoids the main course of decubitus, reduces the rate of postoperative health risk, and speeds up the patient recovery.



Patient would feel creditable comfortable thanks to extra 75mm thick padding of a double-layer construction inside.



The outstanding GORE-TEX material has the following characters: breathable, waterproof, fireproof, anti-static and match the hygiene standard.



Since normal affixed or sewed cushion can't be directly washed, how to clean the blood and fluid stains which will be a hidden danger for surgeon and patient's health has become a serious problem. Using advanced ultrasonically-welded technology, our excellent waterproof cushion will never have such worries.



Foot break provides maximum safety and easy use.



Battery power can support 50-80 operations. Backup control keys in case of hand control failure. Optional stainless steel cover.


Backup control
Longitudinal shift*
up/down
Trendelenburg
Reverse Trendelenburg
Turn left/right
Back plate up/down

*Optional function

HyBase 3000

Electro-Hydraulic Operating Table


Technical Specifications

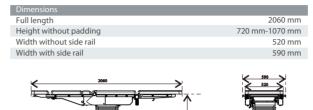


	728-1070	
Frendelenburg position		
Trendelenburg		Max. 30°
Reverse trendelenburg		Max. 30°




Posición cóncava	
Leg plates up	Max. 20°
Leg plates down	Max. 90°
	# K
	20°

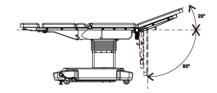


Reflex Position	
Reflex Position	110°
Trendelenburg	25°
Mounting point up	45°

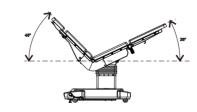

Longitudinal slide(optional)	
Longitudinal slide	Max. 300 mm

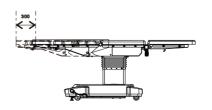

Lateral tilt			
Left tilt			Max. 25°
Right tilt			Max. 25°
	·**	1	

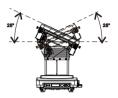

Back plate Back plate up Back plate down	Max. 80° Max. 40°
J	
**	

Head plate	
Head plate normal travel	Max. 45°
Head plate extended travel	Max. 90°
45'	

Flex Position	
Flex Position	220°
Reverse Trendelenburg	20°
Mounting point down	20°



Body elevator(optional)	
Body elevator up	Max. 120 mm

